


*Innovative Patient Care Products  
Designed for Advanced Care*



**DM**  
SYSTEMS  
I N C



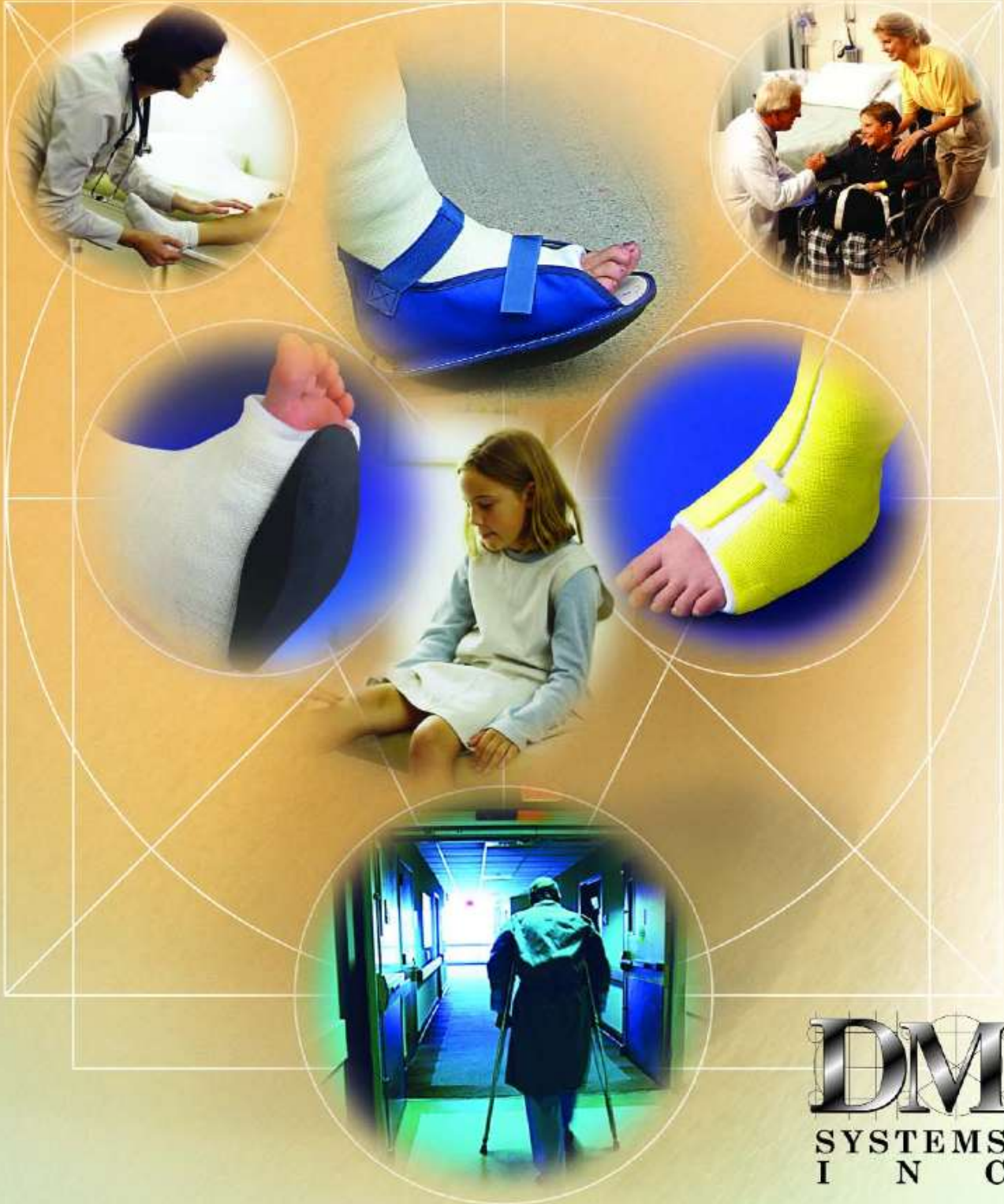

## GaitKeeper™ CAST SHOE

### Removable air-cushioned cast protection

GaitKeeper Cast Shoe with cushioned air sole design prevents excessive walking cast wear and reduces impact loading of the injured leg. An enhanced rocker bottom permits a more natural gait while the flexible sole, designed with 45 pockets of air, delivers superior comfort for the cast wearer.



*Cutaway diagram illustrates how air sole design reduces impact.*

Patent Pending  
Suggested Code:  
L3260



- Easy-on/Easy-off for the convenience of the cast wearer
- Held firmly in place by elastic straps across the top
- Ideal for long or short casts
- Available in 5 sizes—Extra Small, Small, Medium, Large, Extra Large
- Limited three-month/single patient warranty
- Latex-free



## CastWalker® CAST SOLE

### Fixed cast protection designed for improved gait

Increased patient mobility and comfort are benefits of CastWalker Cast Sole, an innovative, effective option for use on walking casts. CastWalker allows for a more normal heel-toe motion so wearing a cast becomes less of an impediment to an active lifestyle.



Patent No. 5,433,695  
Suggested Code: L3260



- Applies easily to the bottom of a walking cast with fast drying adhesive (included)
- Universal design for use on either right or left foot
- Available in 3 sizes—Large, Medium, Small
- Latex-free



# Two Ways To Ease The Discomfort Of A Walking Cast



# CastWedge™

## CAST ADJUSTER

### Accommodates swelling and aids in realignment

To prevent the discomfort of post-trauma swelling, CastWedge permits uni- or bi-valving the cast after the initial injury or surgery. Maintains a precise opening in the cast. For realignment, use CastWedge to make exact adjustments to casted fractures without the inconvenience or cost of replacing the original cast.



- Made of high-impact, non-porous, radiolucent material that will not interfere with X-rays
- Lightweight, inexpensive and available in 5 sizes—10mm, 15mm, 20mm, 25mm, and 30mm
- Use with plaster or fiberglass casts
- Starter set of 20 wedges each of five sizes (100 total) with replacement packages of individual sizes sold separately



Patent No.  
4,372,300

Suggested Code:  
Q4050C



**DM Systems, Inc.**  
1316 Sherman Avenue  
Evanston, IL 60201 USA  
800.254.5438 • Fax: 847.328.9561  
E-mail: [info@dmsystems.com](mailto:info@dmsystems.com)  
[www.dmsystems.com](http://www.dmsystems.com)